



## The Most Efficient Enterokinase for Fusion Tag Cleavage

If you need to cleave your fusion protein with enterokinase, here's good news. Invitrogen's recombinant EnterokinaseMax™ (EKMax™) provides you with unsurpassed digestion efficiency. With EKMax™, you can completely cleave fusion tags using minimal amounts of enterokinase. In addition, the highly specific EK-Away™ resin effectively removes enterokinase from your digestion mixture. With EKMax™ and EK-Away™ you can obtain pure native protein for downstream experiments.

### Why EKMax™?


EKMax™ is a highly purified recombinant enterokinase. It recognizes the specific amino acid sequence Asp-Asp-Asp-Asp-Lys and cleaves following the lysine residue. EKMax™ is the most efficient enterokinase available, because it offers exceptional activity and specificity.

**EKMax™ is the most active enterokinase around.** Almost all commercially available enterokinases are porcine or bovine holoenzymes. EKMax™ is a recombinant preparation of the active catalytic subunit of bovine enterokinase (1). This means that with EKMax™ you use less enzyme than other commercial preparations and still achieve maximum digestion efficiencies.

**EKMax™ is ultra specific.** EKMax™ is produced in the yeast *Pichia pastoris*. This unique expression system secretes high levels of correctly processed enterokinase into the growth medium so EKMax™ contains almost no foreign proteins right from the start. This means EKMax™ has a higher specific activity and there is no unwanted cleavage caused by contaminants.

**EKMax™ really works.** To demonstrate how well EKMax™ works, we've used EKMax™ to digest a thioredoxin-chloramphenicol acetyl transferase (CAT) fusion protein. Figure 1 shows that protein is completely cleaved into thioredoxin and CAT with just 0.6 units\* of EKMax™.

Figure 1 - Digestion of thioredoxin-CAT fusion protein with EKMax™

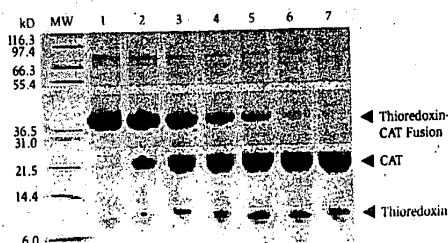

Each reaction contains 10 µg of partially purified thioredoxin-CAT fusion protein and varying amounts of EKMax™. Reactions were incubated at 37°C for 16 hours and analyzed on a Coomassie® blue-stained 10-20% Thine-SDS gel. Units of EKMax™ used per reaction are listed below.

Lane 1: No EKMax™	Lane 5: 0.2 units
Lane 2: 0.05 units*	Lane 6: 0.4 units
Lane 3: 0.1 units	Lane 7: 0.6 units
Lane 4: 0.15 units	


\* One unit of EKMax™ is the amount of enzyme required to digest 20 µg of a thioredoxin-CAT fusion protein to 90% completion in 16 hours at 37°C. One EKMax™ unit is equivalent to ~190 trypsinogen activation units.

## EK-Away™ for rapid removal of EKMax™

If removal of EKMax™ from your native protein is critical for subsequent analysis, try EK-Away™. EK-Away™ is an agarose based resin that provides a quick and efficient way to separate enterokinase (EKMax™ or enterokinase holoenzymes)

from a mixture of proteins by binding to the enzyme's catalytic site. It takes only 5 minutes for the EK-Away™-enterokinase complex to form. Following brief centrifugation, greater than 99% of the enterokinase can be removed (figure 2).

Figure 2 - Western blot of EKMax™ bound to EK-Away™

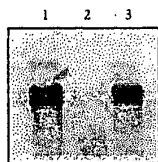

One hundred and fifty microliters of EK-Away™ was added to 500µl cleavage mixture containing 5 units of EKMax™. The EK-Away™ resin from the binding reaction was spun down, suspended in Laemmli loading buffer and boiled for 1 minute to release EKMax™. Ten microliters of the boiled supernatant, along with 10 µl from the supernatant of the binding reaction, were separated by SDS-PAGE. The gel was blotted and probed using an anti-enterokinase antibody.

Lane 1: 0.1 units of EKMax™ (control)  
Lane 2: 10 µl of binding reaction supernatant  
Lane 3: 10 µl of supernatant from boiled resin

## Quality guaranteed

To guarantee consistent quality results with EKMax™ and EK-Away, each product is stringently tested. Our quality control standards require that 1 unit of EKMax™ digests 20 µg of a thioredoxin-CAT fusion protein to 90% completion in 16 hours at 37°C. Each lot of EKMax™ is also incubated

with azocasein to ensure that there is no non-specific protease activity. Each lot of EK-Away™ resin is tested by a sensitive fluorometric assay to ensure that less than 1% of enterokinase activity remains after the binding reaction.

## Ordering information

EKMax™ and EK-Away™ together represent the most efficient system available for enterokinase cleavage and removal. EKMax™ is supplied in 250 and 1,000 unit sizes and comes with 10X reaction buffer and a comprehensive user manual.

EK-Away™ is offered in two convenient sizes that include the EK-Away™ resin, 10X binding buffer, and 10X stripping buffer. To get your ultra pure recombinant native protein, call Invitrogen and order EKMax™ and EK-Away™ today!

Description	Quantity	Cat. no.
EnterokinaseMax™ (EKMax™)	250 units	E180-01
	1,000 units	E180-02
EK-Away™	7.5 ml*	R180-01
	30 ml*	R180-02

\* 7.5 ml EK-Away™ resin will remove 250 units of EKMax™; 30 ml will remove 1,000 units of EKMax™.

### Reference:

1. La Vallie, E.R. *et al.* (1993) *J. Biol. Chem.* 268: 23311-23317.



B-01E030231 020701

1600 Faraday Avenue • Carlsbad • CA 92008  
Phone: 760 603 7200 Fax: 760 602 6500  
Toll Free: 800 955 6288 • [www.invitrogen.com](http://www.invitrogen.com)

This product is sold under patent license from Genetics Institute, Inc. License for commercial use may be obtained directly from Genetics Institute. All prices are in U.S. dollars and are subject to change. For research use only.  
© Invitrogen Corporation 2001. Reproduction forbidden without permission.